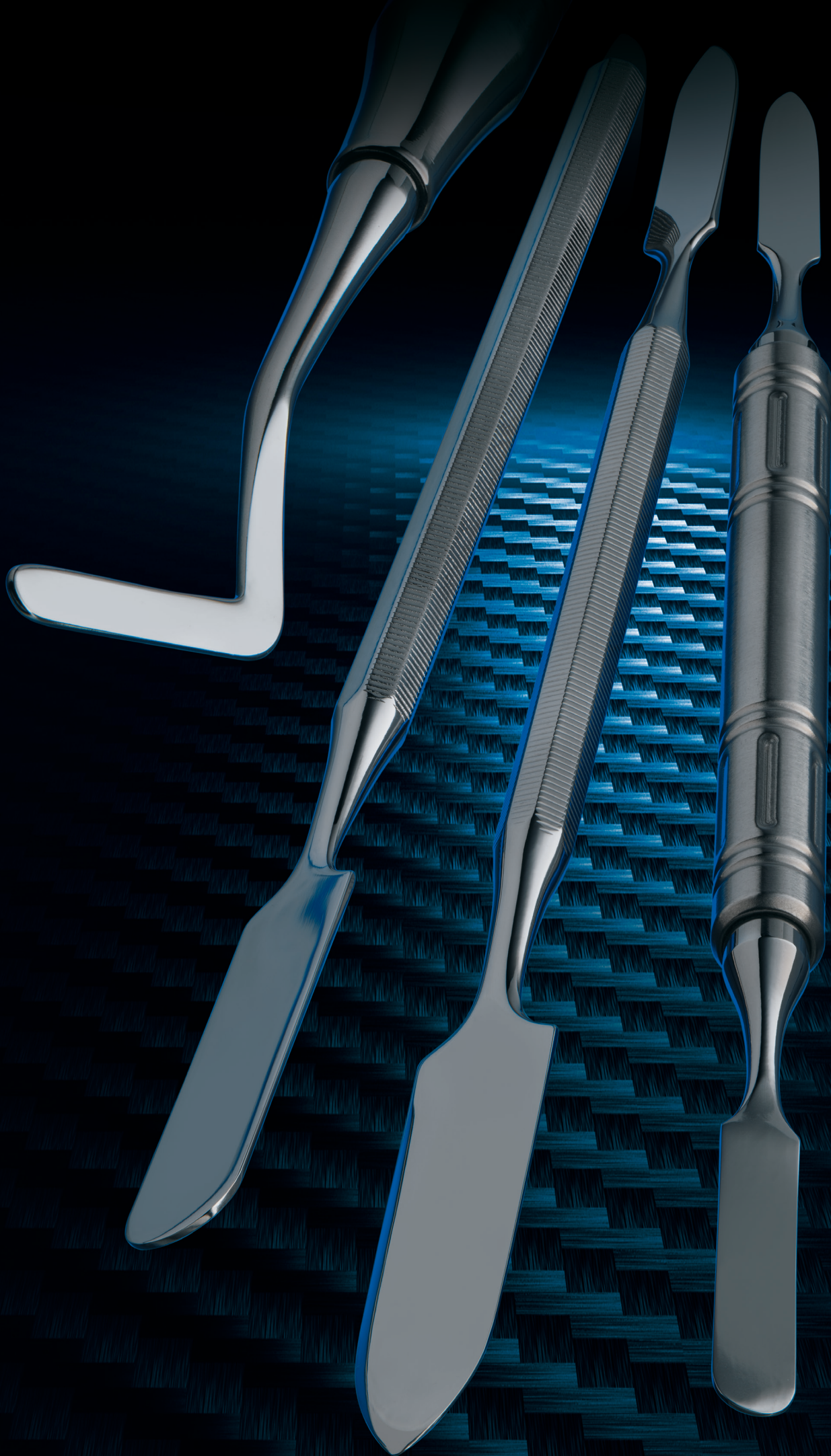
03

05

Spatole per cemento e Heidemann

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Cement spatulas and Heidemann

Heidemann spatulas and cement spatulas are dental instruments used to mix and apply dental cements, materials used to attach reconstructions such as crowns, bridges, and infills on teeth, to fill cavities, to cement orthodontic brackets and bands, and to perform other dental procedures. These instruments have flat stainless-steel blades designed to be smooth and nonporous to prevent dental cement from sticking to the blade. Once the cement has been mixed properly, they are used to apply the material to the tooth or onto the reconstruction for placement. The flat blade of the spatula enables accurate cement placement and helps ensure that it gets evenly and smoothly applied.

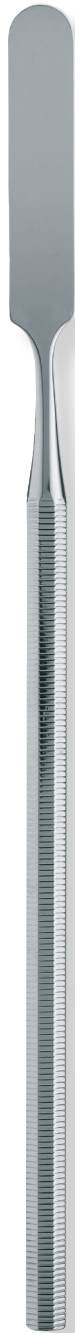
fig. **2**
0300.E0.02

fig. **0A**
0303.E0.0A

fig. **1**
0303.E0.1

fig. **2**
0303.E0.2

fig. **3**
0303.E0.3

fig. **4**
0303.E0.4

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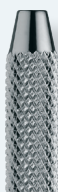
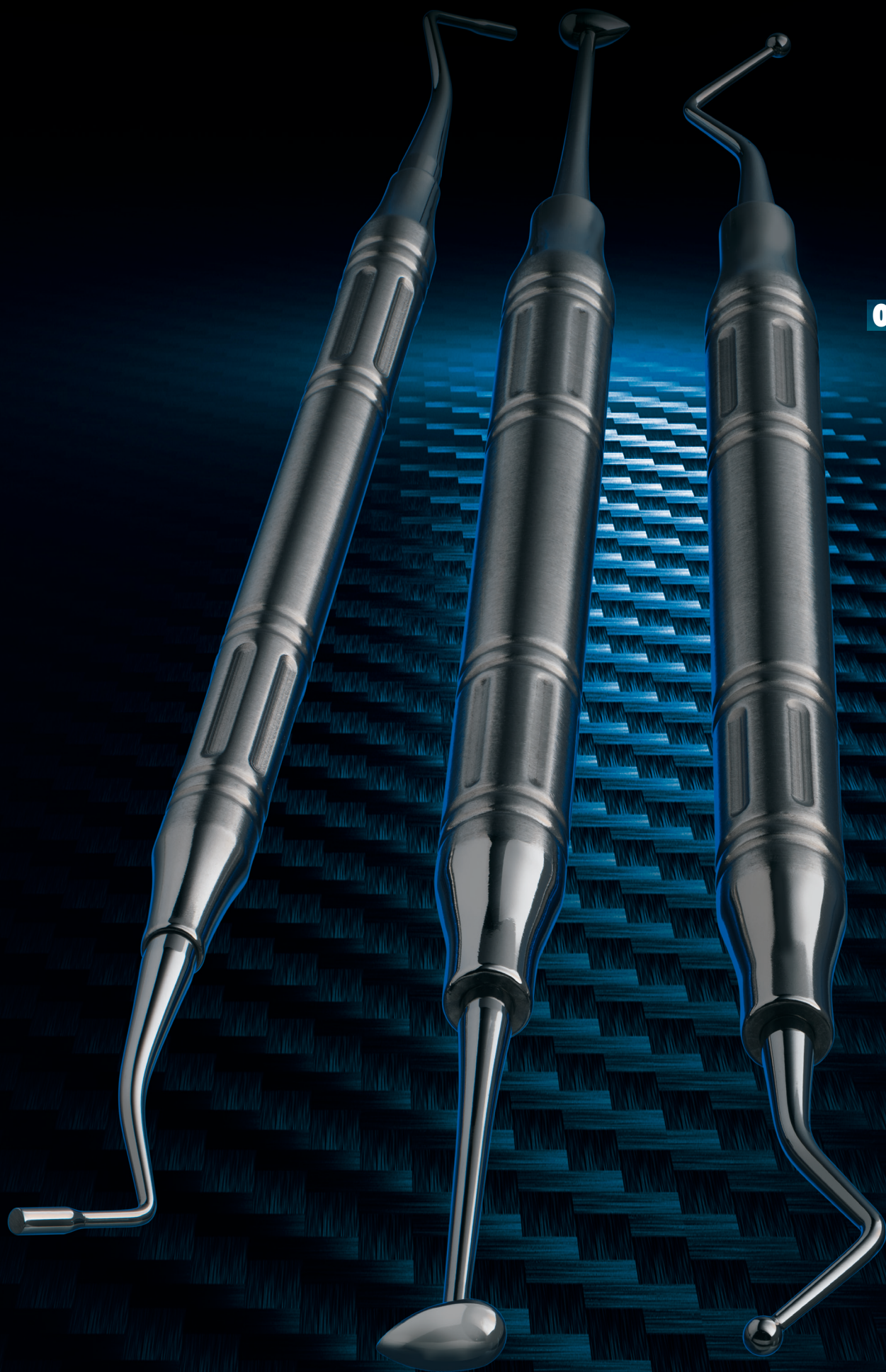
The dental filling procedure is included in a category of practices known as conservative dentistry. This particular branch of dentistry focuses on the restoration of decayed or damaged teeth, once the underlying process of decay has been eliminated. The term "conservative" is used because the tooth is preserved in its original location and undergoes repairs to the greatest extent possible, aiming to restore its structural integrity and prevent extraction. Dental fillings, also referred to as dental restorations, involve the placement of restorative material within a cavity or lesion in order to repair tooth decay.

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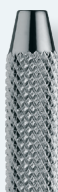
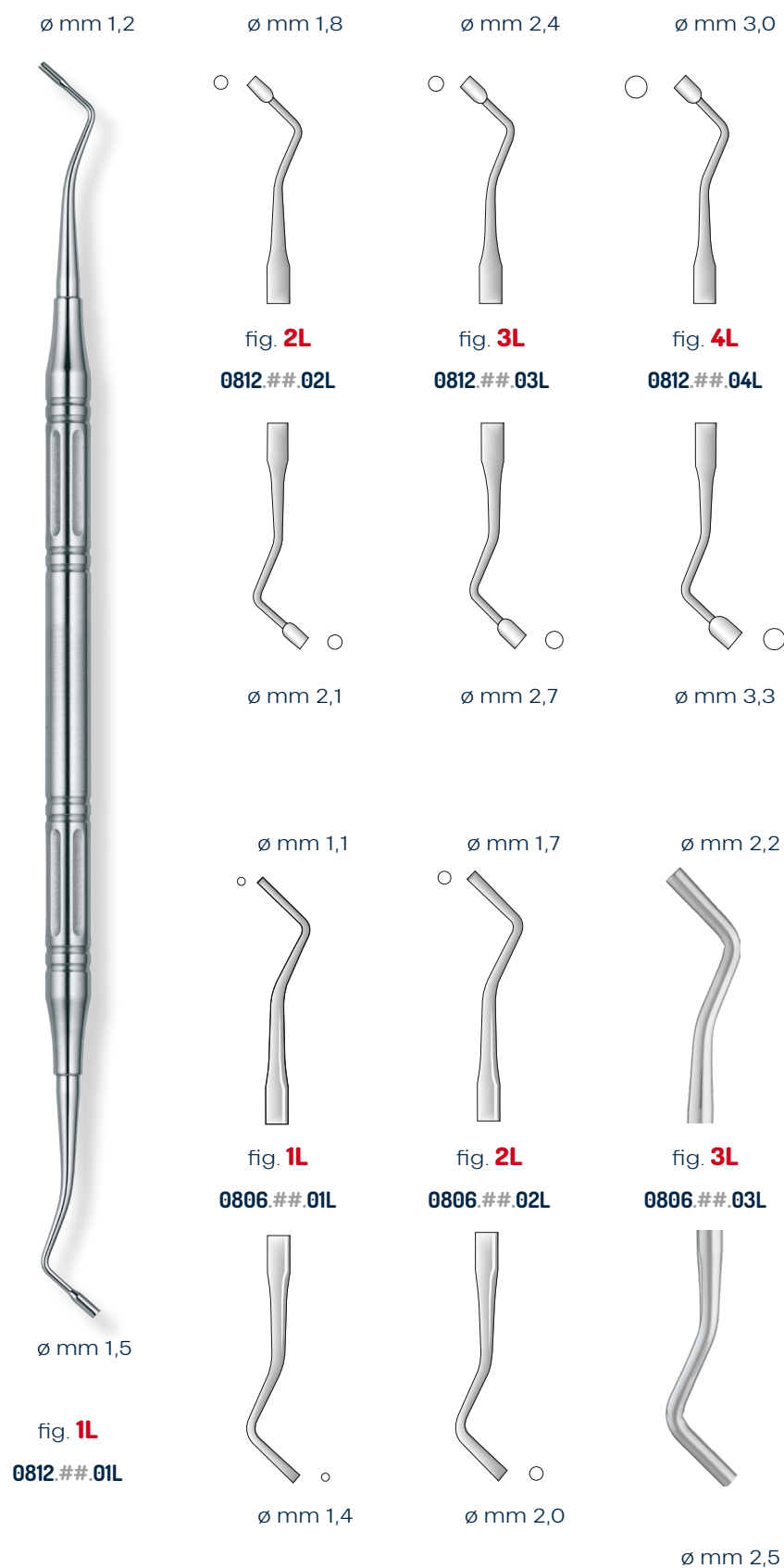
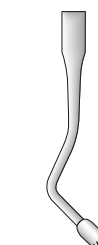
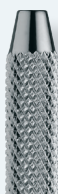
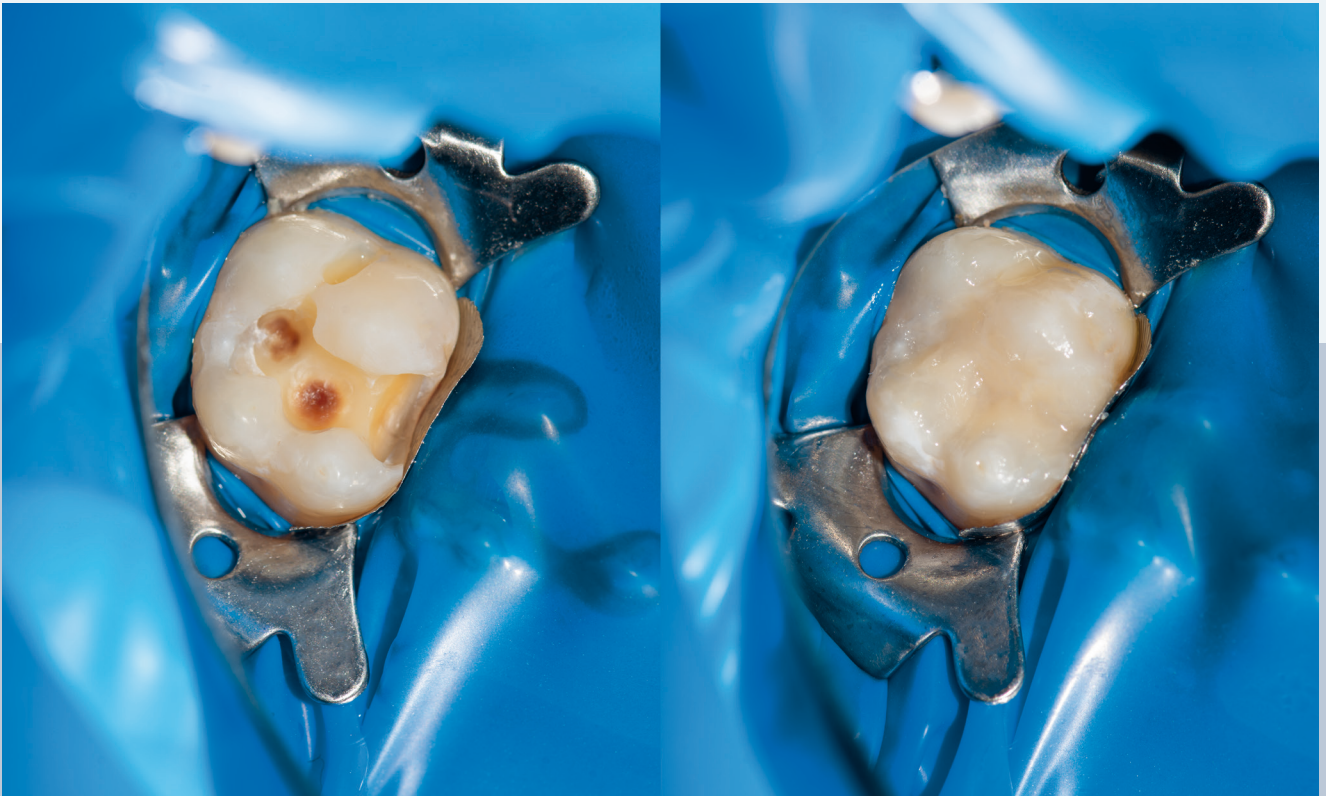
Dental filling instruments are utilized to compress and pack dental filling materials, such as amalgam or composite resin, into a pre-treated and prepared tooth cavity or other dental structure. The primary objective of using these instruments is to ensure the accurate placement, fitting, and packing of the filling material within the cavity or structure. This is crucial in order to prevent the formation of gaps or hollow areas that may result in tooth sensitivity, marginal leakage, or other potential issues.


Burnishers are commonly employed during the final stages of dental restoration procedures or surgeries. Their purpose is to accentuate grooves and eliminate excess fins, while also refining the rough edges following the dental restoration. Additionally, they serve to polish the dental restoration, ensuring a smooth finish, and eliminate any scratches that may have arisen on the amalgam surface during the restoration process.

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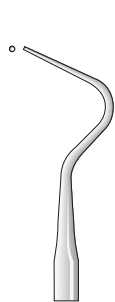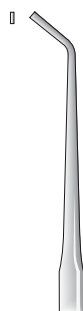
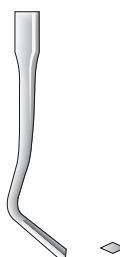
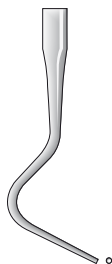
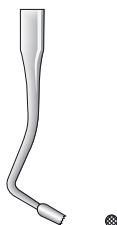
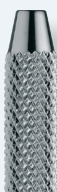
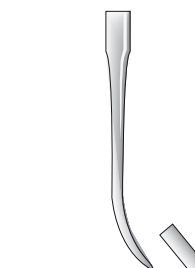
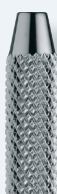
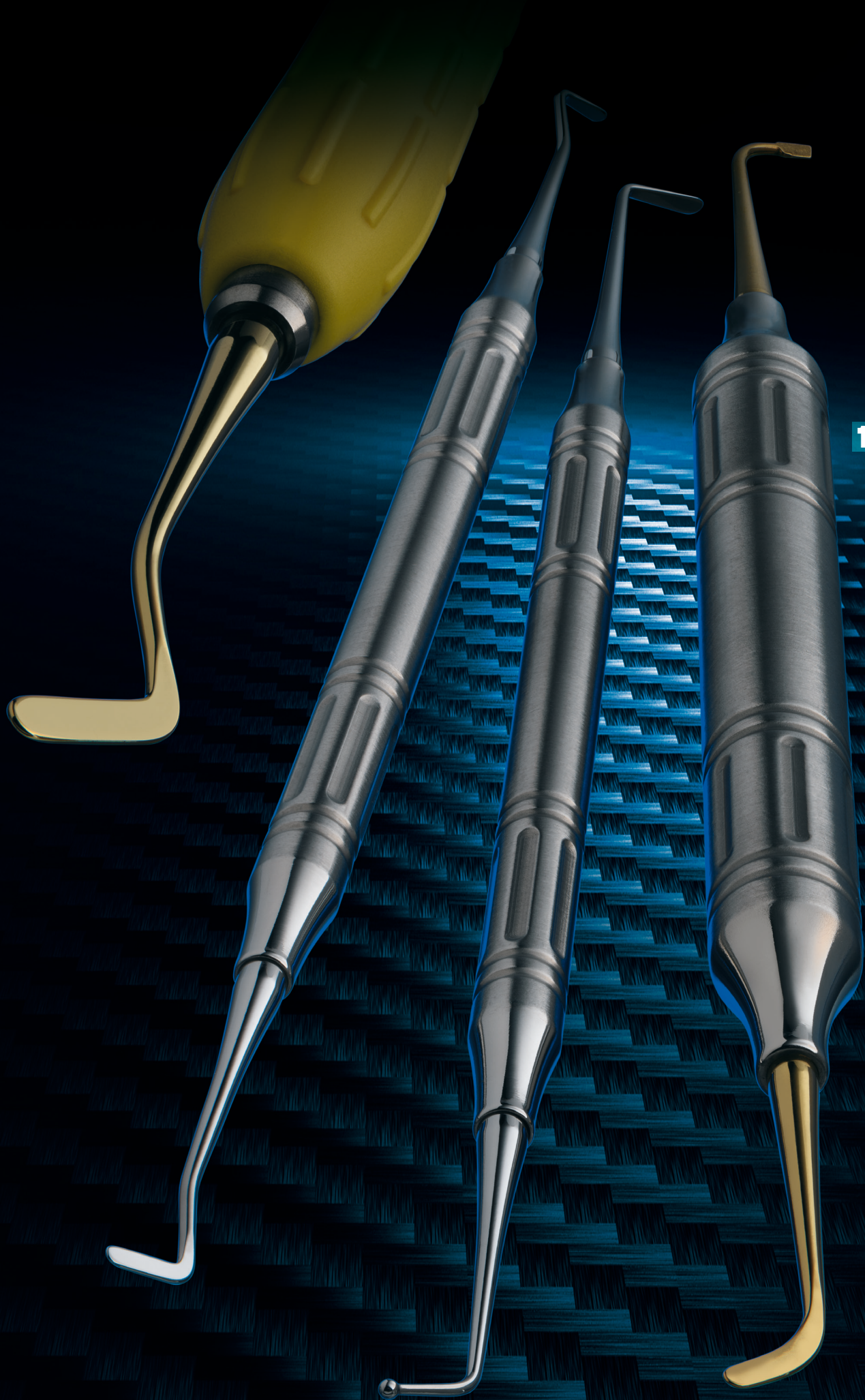
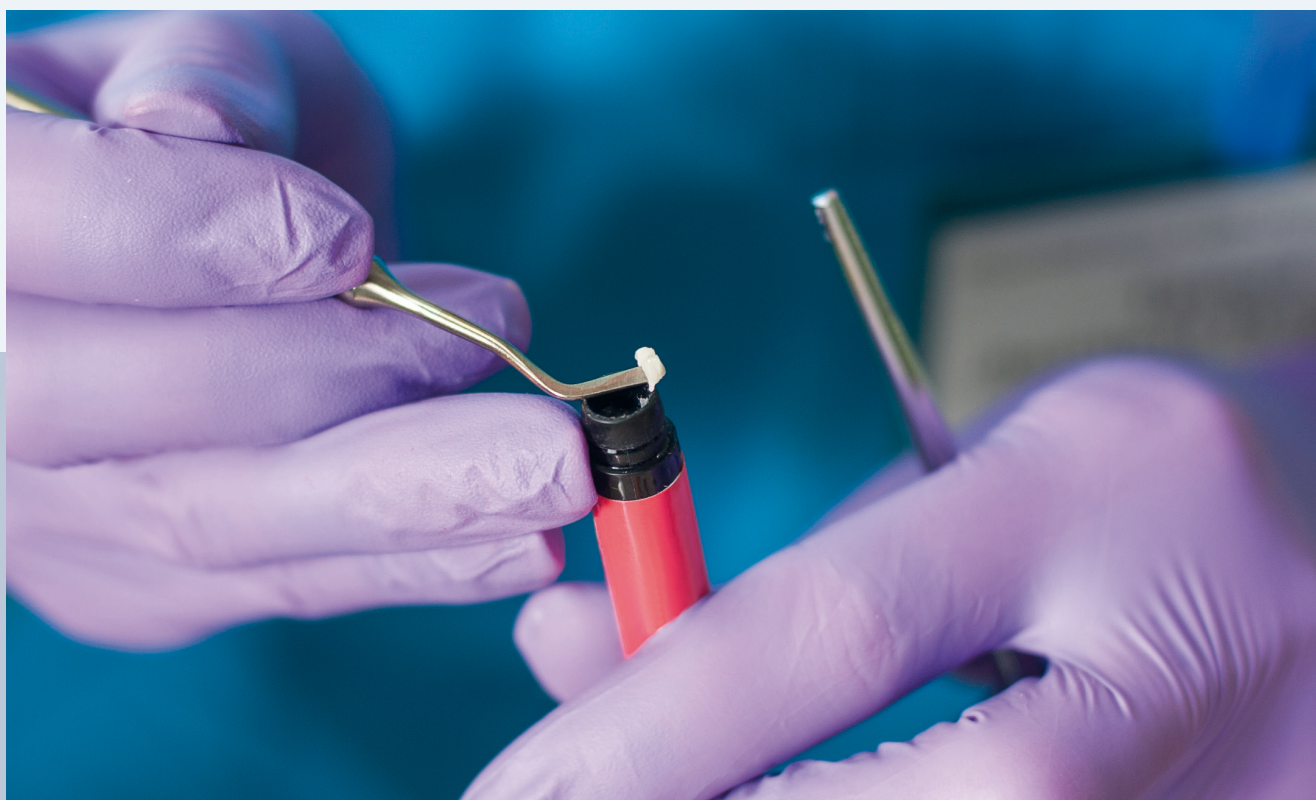
Composite Instruments

Dental composite is a restorative material used in dentistry to repair carious cavities or dental lesions. It is a mixture of acrylic resin and filler particles, such as silica, quartz or glass, and is used to mold and restore the shape and function of the damaged tooth.

Dental composite has several advantages over other restorative materials, such as dental amalgam. First, the color of the composite can be matched to the one of what is left of the tooth, making the filling almost invisible. Second, the composite adheres firmly to the tooth, reducing the risk for the fixture to fracture or fall out.

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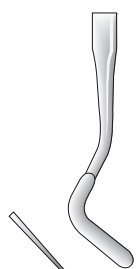
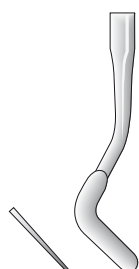
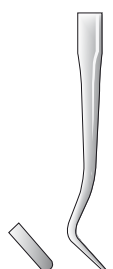
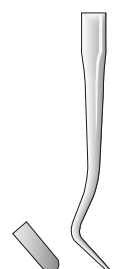
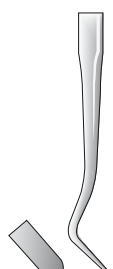
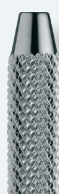
Composite instruments essentially serve three functions:

- 1.** Composite placement - used to place and shape the composite resin material in the pretreated cavity.
- 2.** Mold/shape the composite - used to contour and shape the composite resin after it's been placed.
- 3.** Polish the composite - used to smooth and polish the surface of the composite resin restoration.

fig. **03**
1115.##.03

fig. **6**
1115.##.06

fig. **11**
1115.##.11

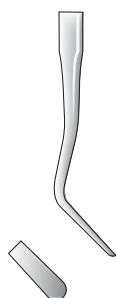
fig. **179**
1115.##.179

fig. **00**
1107.##.00

fig. **33**
1107.##.33

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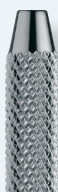
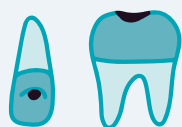
CLASSE V

CLASSE VI


Black's classification of dental caries is one of the most widely used in dentistry. It was developed by American dentist G.V. Black in 1908 and still remains a benchmark for the diagnosis and treatment of dental caries.

Black's classification divides carious cavities into five to six classes, based on the anatomical location on the tooth and their extent:

- **CLASS I:** caries affects the occlusal surface at the level of pits and fissures, the buccal and lingual/palatal surfaces at the level of molars and premolars, or the lingual pits of incisors and canines.
- **CLASS II:** Caries is located in-between teeth (interproximal walls) of molars and premolars. In most cases, class II carious lesions also affect the occlusal surface.
- **CLASS III:** Caries affects the interproximal walls of incisors and canines without involvement of the incisal edge (corner of the tooth).
- **CLASS IV:** caries involves incisors and canines simultaneously affecting their interproximal walls and incisal edge
- **CLASS V:** caries involves the collars (area of the tooth close to the gum) of both anterior and posterior teeth
- **CLASS VI:** Caries affects the incisal margin of the anterior elements or the cusp tips of canines and molars and premolars.

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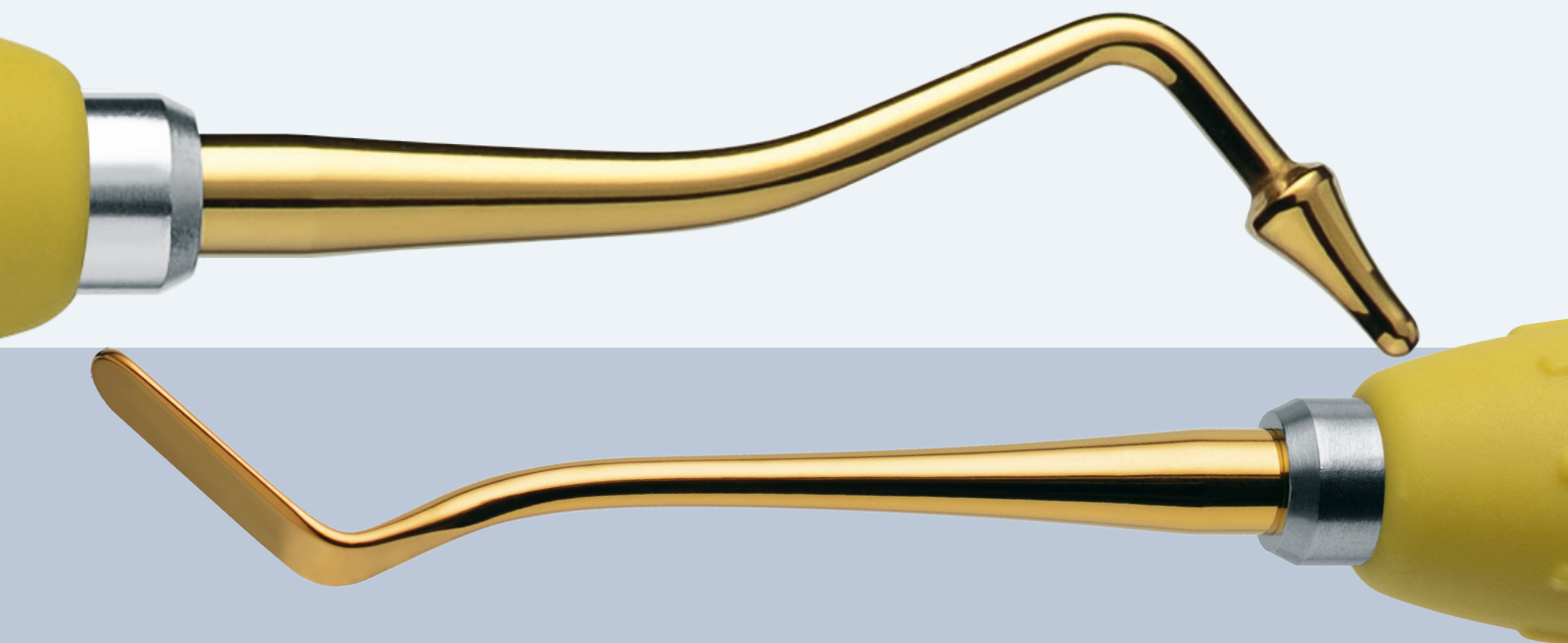Zirconium nitride (ZrN) is a hard, wear-resistant, biocompatible material that is used to improve the performance and durability of dental instruments.

The coating is created using a process called physical vapor deposition (PVD), which involves the deposition of zirconium nitride atoms on the surface of the instrument. The resulting coating is a thin, smooth, non-stick layer that improves surface hardness and reduces friction, making the instrument more resistant to wear and corrosion.

Some advantages of using ZrN coating on dental instruments are:

- 1.** Increased durability and longevity of the instrument.
- 2.** Reduced friction and improved cutting performance.
- 3.** Significantly improved surface non-stickiness resulting in better workability of the composite.
- 4.** Improved wear and corrosion resistance.
- 5.** Improved biocompatibility and reduced risk of allergic reactions.

Strumenti per composito rivestiti in Nitruro di Zirconio - ZrN Zirconium Nitride Coated Composite Instruments - ZrN

11
59

fig. 1

1110.##.01TN

fig. 2

1110.##.02TN

fig. 3

1110.##.03TN


fig. 4

1110.##.04TN

fig. 10

1111.##.10

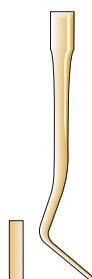Goldstein fig. 1

1111.##.01


Goldstein fig. 2

1111.##.02

Goldstein fig. 3

1111.##.03

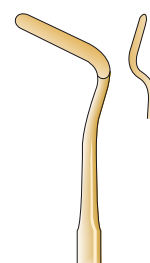
Goldstein fig. 4

1111.##.04

Goldstein fig. 5

1111.##.05
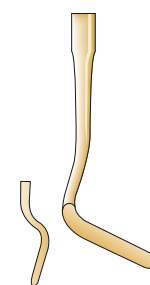
##

A0

G0

DT06

DT08

DT10

EVO MY

Manici disponibili

Available handles


11
60



fig. 10
1110.##.10TN

Felt fig. 1
1110.##.FELT1

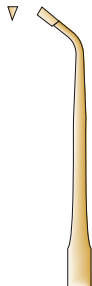Felt fig. 2
1110.##.FELT2

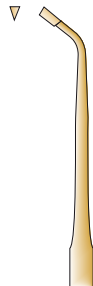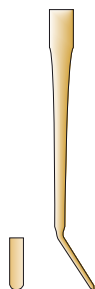Felt fig. 3
1110.##.FELT3

Felt fig. 4
1110.##.FELT4

Felt fig. 5
1110.##.FELT5

Felt fig. 6
1110.##.FELT6

Felt fig. 7
1110.##.FELT7

##

A0

G0

DT06

DT08

DT10

EVO MY

**Manici disponibili
Available handles**

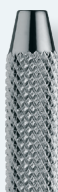
fig. **1 Applicator**
118.##.01

fig. **2 Condenser**
118.##.02

fig. **3 Modelling**
118.##.03

fig. **4 Control**
118.##.04

fig. **5 Remover**
118.##.05