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**Periodontal curettage or scaling** is a dental treatment that involves cleaning the roots of the teeth to remove tartar, bacterial plaque, and infected or inflamed tissues from the area around the roots.

During the periodontal curettage process, the dentist or dental hygienist uses special instruments like curettes and scalers to scrape away the tartar and plaque that have accumulated below the gum line and along the roots of the teeth. This helps to remove sources of inflammation and infection, promoting the healing of the gums and improving periodontal health.

**Curettes** are manual instruments designed to remove tartar and plaque from the root surface of the teeth and periodontal pockets.

They have rounded tips and a curved cutting edge. This shape allows for safe access to deep periodontal pockets without damaging the soft tissues.

**Scalers** are manual instruments primarily used for the removal of supragingival tartar.

They have pointed tips and sharp edges on both sides of the blade. This design is ideal for scraping away tartar above the gum line.

# Gracey Curettes

**34**

137

fig. **3/4**

3400.##.03/04

fig. **5/6**

3400.##.05/06

fig. **7/8**

3400.##.07/08

fig. **9/10**

3400.##.09/10

fig. **11/12**

3400.##.11/12

fig. **13/14**

3400.##.13/14

fig. **15/16**

3400.##.15/16

fig. **17/18**

3400.##.17/18

**Standard  
Blade**

**Standard  
Shank**

**Standard**

fig. **1/2**

3400.##.01/02

Gracey curettes disponibili in versione rigida con codice **3402**  
Gracey curettes with rigid tips availables with code **3402**

**Manici disponibili**  
**Available handles**

##

A0

G0


P0

DT06

DT08

DT10

EVO - M

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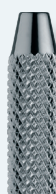
Acciaio INOX • Titanio • DLC

**Available tips in**

Stainless Steel • Titanium • DLC coating

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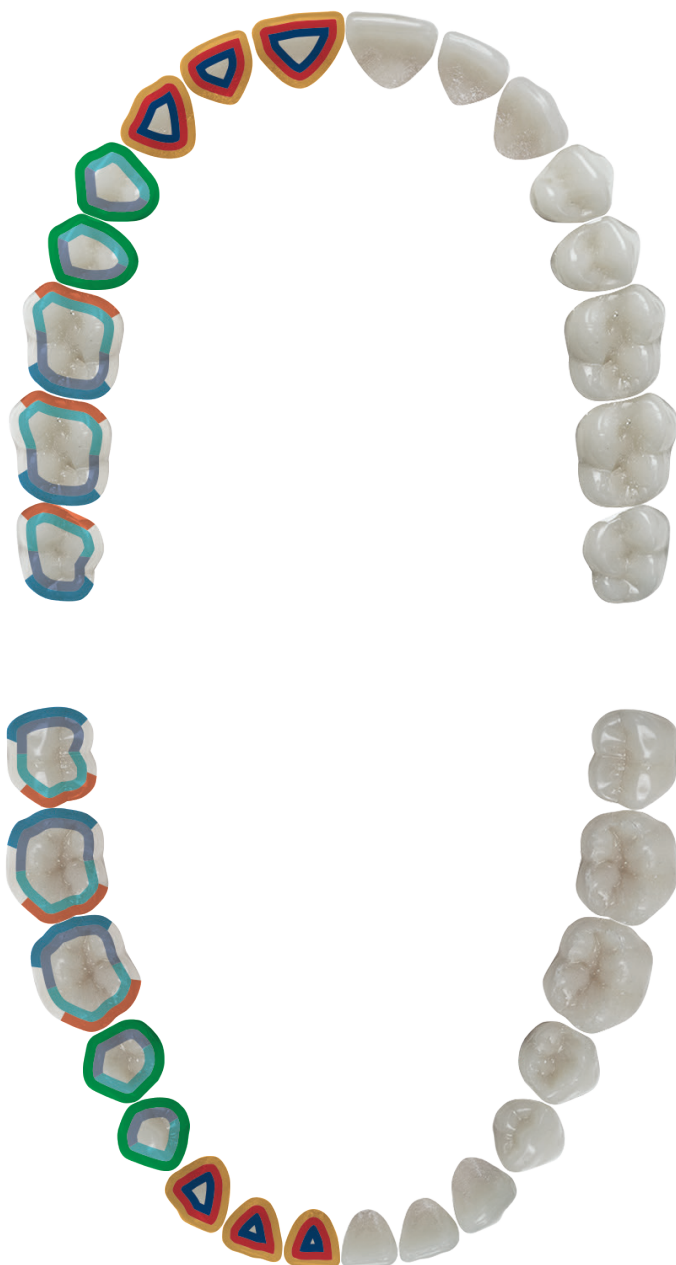
## Area-specific Curettes

Area-specific curettes are designed with blades of different shapes and angles depending on the teeth or distinct areas within the mouth where they will be used. They are customized to fit the anatomy of specific teeth, making them more efficient for scaling and root planing in those areas. They offer better adaptability and efficiency in designated areas due to their specific design. The most important and well-known among these are certainly the **Gracey**.

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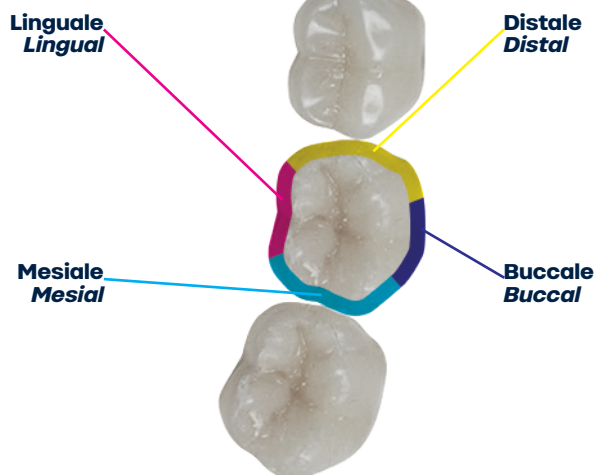**The blades of Gracey curettes** are angled at 60 or 70 degrees relative to the lower shank, and this specific inclination allows for better adaptability to the tooth surface and enables more precise removal of tartar below the gum line.

**Each blade of the Gracey curettes** has only one cutting edge, which is the lower edge when the terminal shank is held vertically. This helps reduce the risk of damaging soft tissues and facilitates control of the instrument during procedures.

**The blades of Gracey curettes** are curved in two planes: upward and sideways, to allow better adaptability to the convex and concave surfaces of the tooth root. They are also particularly effective in cleaning and treating periodontal pockets and are ideal for root planing.

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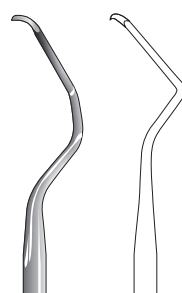
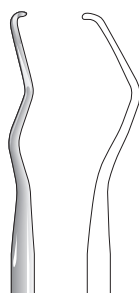
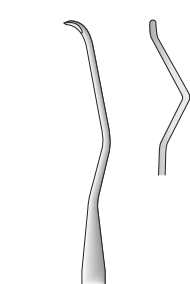
The **Short-Type** version is designed to facilitate insertion into shallow pockets, areas with limited access such as narrow interdental spaces, and to reach smaller roots, furcations, and grooves. Compared to the standard form, the terminal shank is 3 mm longer, and the blade is thinner and 50% shorter.

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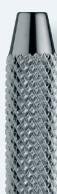
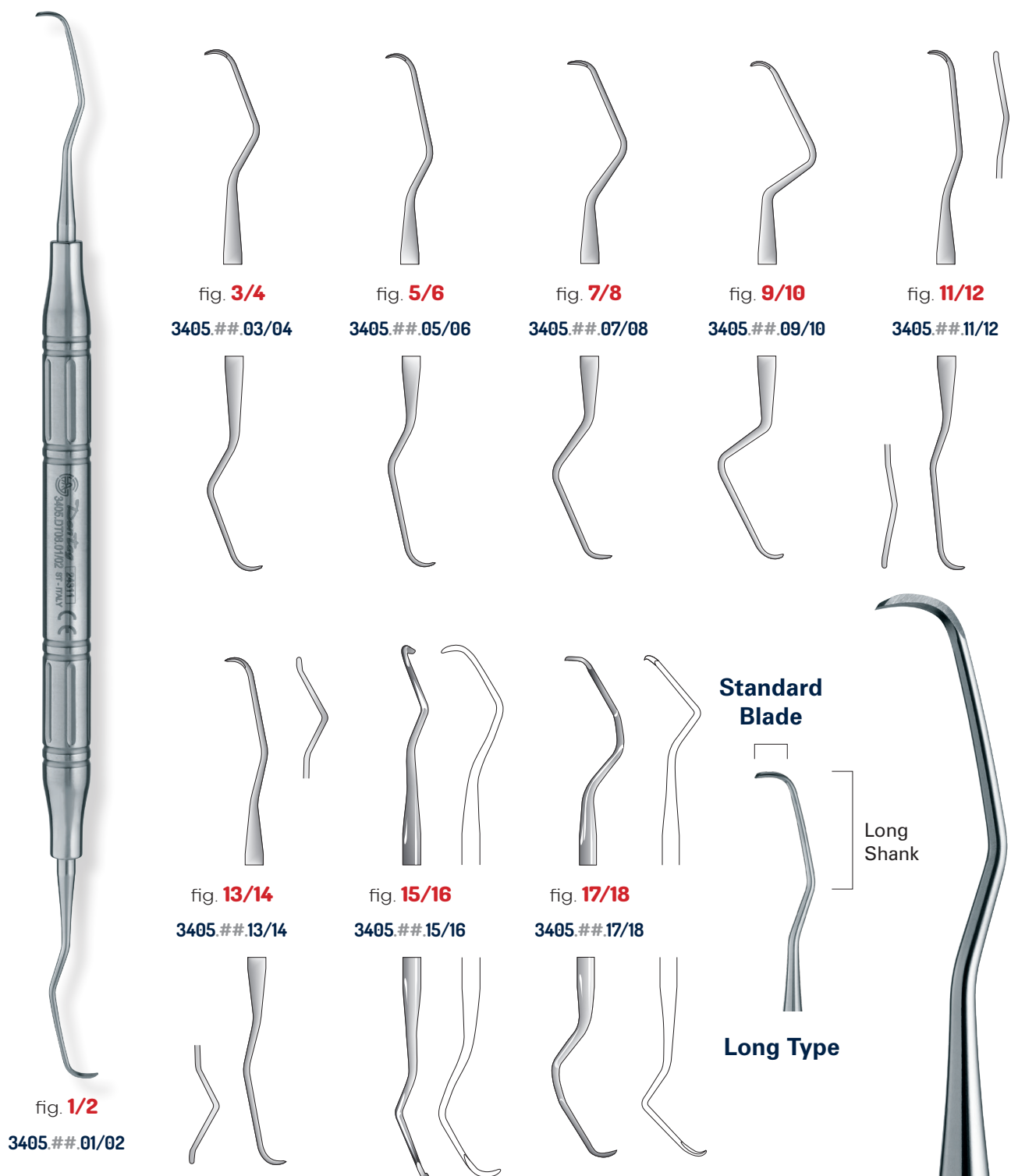
The **Long-Type** version is designed to facilitate subgingival insertion, especially for deeper periodontal pockets and furcations. By reducing tissue distension, they aid in the removal of deposits around the crown and allow better access to the roots. Compared to the standard form, **the terminal shank is 3 mm longer, and the blade is 10% thinner**.

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## **What's EVO handle**

*We are always very attentive to changes and trends in a market that evolves rapidly and sometimes abruptly. The results of this analysis led to the design and production of a handle with several new features and improvements compared to other similar products on the market.*

### **Material**

*In most cases, similar products have silicone handles which, while being an excellent material resistant to high temperatures, have tactile characteristics that, in the case of prolonged use, lead to physical tension in the carpal tunnel resulting in fatigue and pain. For this reason, we chose a different path, using a **very lightweight composite material** (11 g), **resistant to stress**, with a "dry" physical contact, non-fatiguing, **secure grip**, resistant to high temperatures, and therefore **sterilizable in an autoclave**.*

### **Shape**

***Grip diameter of 10.5 mm**, to minimize carpal tunnel problems caused by prolonged use of instruments with smaller grips*

Made in Italy  
Dentag

EVO

Made in Italy  
Dentag

EVO

fig. **1/2**  
**3400.MY.1/2**

fig. **3/4**  
**3400.MR.3/4**

fig. **5/6**  
**3400.MB.5/6**

fig. **7/8**  
**3400.MG.7/8**

fig. **9/10**  
**3400.MW.9/10**

fig. **11/12**  
**3400.MK.11/12**

fig. **13/14**  
**3400.MV.13/14**

fig. **15/16**  
**3400.MO.15/16**

fig. **17/18**  
**3400.ML.17/18**

# Dentag

# Evo similar in body but completely different in Soul

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
## Why choose EVO?

**State of the Art:** all similar instruments on the market are essentially produced through two molding processes. One involves inserting an aluminum or steel tube into the mold, onto which thermoplastic injection is performed. Subsequently, the tips are inserted into the aluminum tube.

Another system involves inserting a finished instrument with a central core of reduced diameter (2.0 – 3.2 mm) into the mold, followed by thermoplastic injection.

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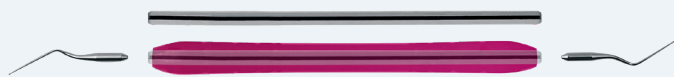## The EVO / Dentag® solution

The way we manufacture the handle is by injecting the composite plastic material in the mould on which two ends we placed two small stainless-steel, shaped bushings. The tips are then inserted in the resulting handle.

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## The difference is evident

**In similar instruments** currently on the market, we find the presence of two distinct materials within the handle (Thermoplastic or Silicone / Aluminum or Steel).

**In our solution,** the inside of the handle is entirely composed of a single composite thermoplastic material. This process eliminates the presence of internal longitudinal metal cores, resulting in significant weight reduction and the elimination of play and cracks due to the contact of very different materials. This also eliminates potential bacterial proliferation problems caused by the creation of air bubbles due to different material expansion gradients during the heating/cooling phases of sterilization cycles.

**Aesthetics:** simple profile, easy to wash and clean. No deep grooves or notches that could cause the accumulation of germs and bacteria. Being a plastic material, it can be colored in various aesthetically pleasing shades, providing advantages for immediate instrument recognition.

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## What's DLC coating

The **DLC** treatment, an acronym for "**Diamond-Like Carbon**", is a coating process that applies a thin layer of carbon with properties similar to those of diamond onto a surface. This coating is known for its excellent mechanical, chemical, and tribological characteristics (related to friction, wear, and lubrication).

## Properties of DLC

**High hardness:** DLC has a hardness that approaches that of diamond, **making surfaces highly resistant to wear**.

**Low friction:** the coating significantly reduces friction between moving surfaces, **improving efficiency and reducing wear**.

**Wear resistance:** thanks to its hardness and low friction, **DLC protects surfaces from abrasion and wear**.

**Chemical resistance:** **DLC is resistant to chemicals and corrosion**, increasing the durability of coated materials.

**Biocompatibility:** **DLC is biocompatible**, meaning **it can be used in medical and food applications without causing adverse reactions**.

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## Why choose DLC coating

**Resharpen dental instruments**, especially smaller ones, **has always been a major challenge** for all users (doctors, dentists, or assistants), because very few can carry out this task, that requires very good manual skills, without compromising the integrity of the instrument itself. Hence, we have developed a new range of instruments that, by using new nano-technologies, through the deposition of extremely hard and resistant ceramic materials, **tackles the root of the problem, and cancels or significantly reduces the need for sharpening**.

To achieve this, we started from raw materials, using only higher quality stainless steels, that underwent rigorous heat treatments and were meticulously sharpened in-house. **Tips were then DLC coated**.

Moreover, the process takes place at low temperatures, exploiting vacuum deposition with the help of Plasma, without affecting the hardening and mechanical characteristics of the steel underneath.

The result is high resistance to abrasion and wear (**High hardness**), low friction coefficient (Non-stick) and excellent chemical inertia (**Corrosion resistance**).

**And, not least, full biocompatibility.**

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
**Universal periodontal curettes** can be used on virtually all tooth surfaces; this is why they are called "Universal." The blade is angled at 90° to the shank, is curved in only one plane, and has both cutting edges. These features make them particularly suitable for the removal of light to moderate amounts of tartar.

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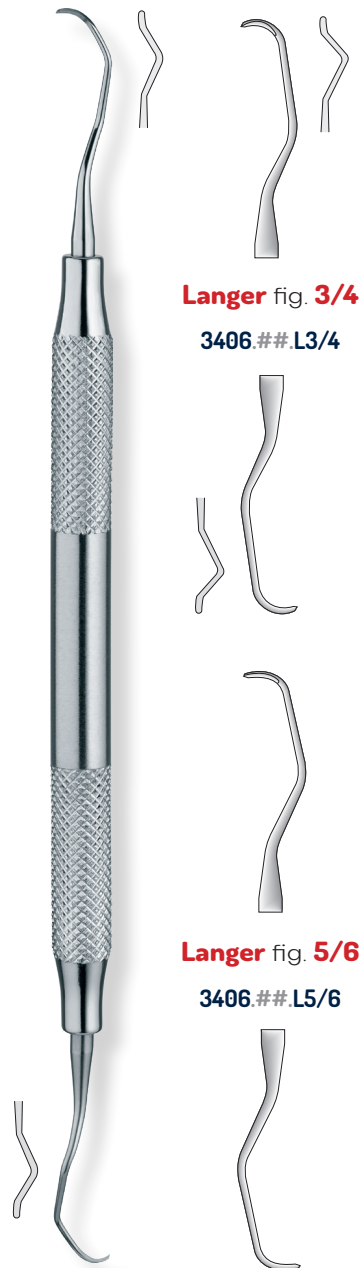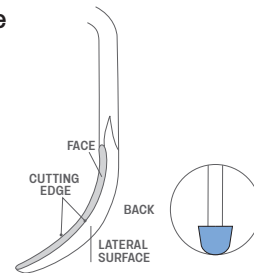
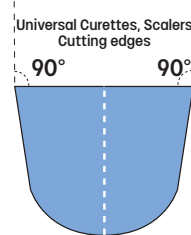
The **Langer curette** combines a universal blade, characterized by two parallel edges and a 90-degree angle between the terminal shank and the blade face, with the shank design of area-specific Gracey curettes. The shapes of the Langer curettes **allow for the removal of tartar and plaque from both mesial and distal surfaces**, simplifying the instrument selection for the hygienist and requiring fewer transitions between instruments. The Gracey-style shank angulation is ideal for hard-to-reach areas, such as deep periodontal pockets and posterior root surfaces. **With only three instruments, it is possible to access the entire mouth.**

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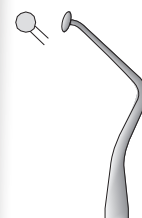
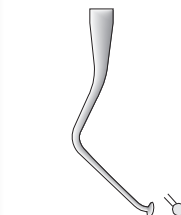
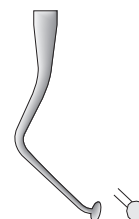
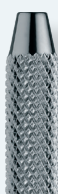
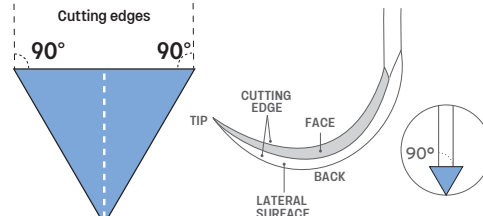
**Scalers are instruments designed for the removal of supra-gingival calculus and surface stains with their sharp and pointed tips.** They are primarily used for the rapid removal of large deposits of calculus visible above the gum line. They are very effective in removing hard and encrusted calculus above the gum line. They are less suitable for subgingival use due to their pointed tip, which could cause trauma to the soft tissues.

**Sickle Scalers:** These have a pointed tip and are mainly used to remove large deposits of calculus from the crowns of teeth, particularly effective in interproximal areas and on supra-gingival surfaces.

**Jacquette Scalers:** These feature a triangular-shaped tip and are designed to access deep periodontal pockets and remove calculus and plaque from root surfaces.

**Hoe Scalers:** These have blades similar to a gardening hoe, with a broad, flat blade, and are used to remove calculus and plaque from the tooth surfaces below the gum line.

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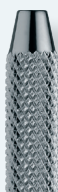
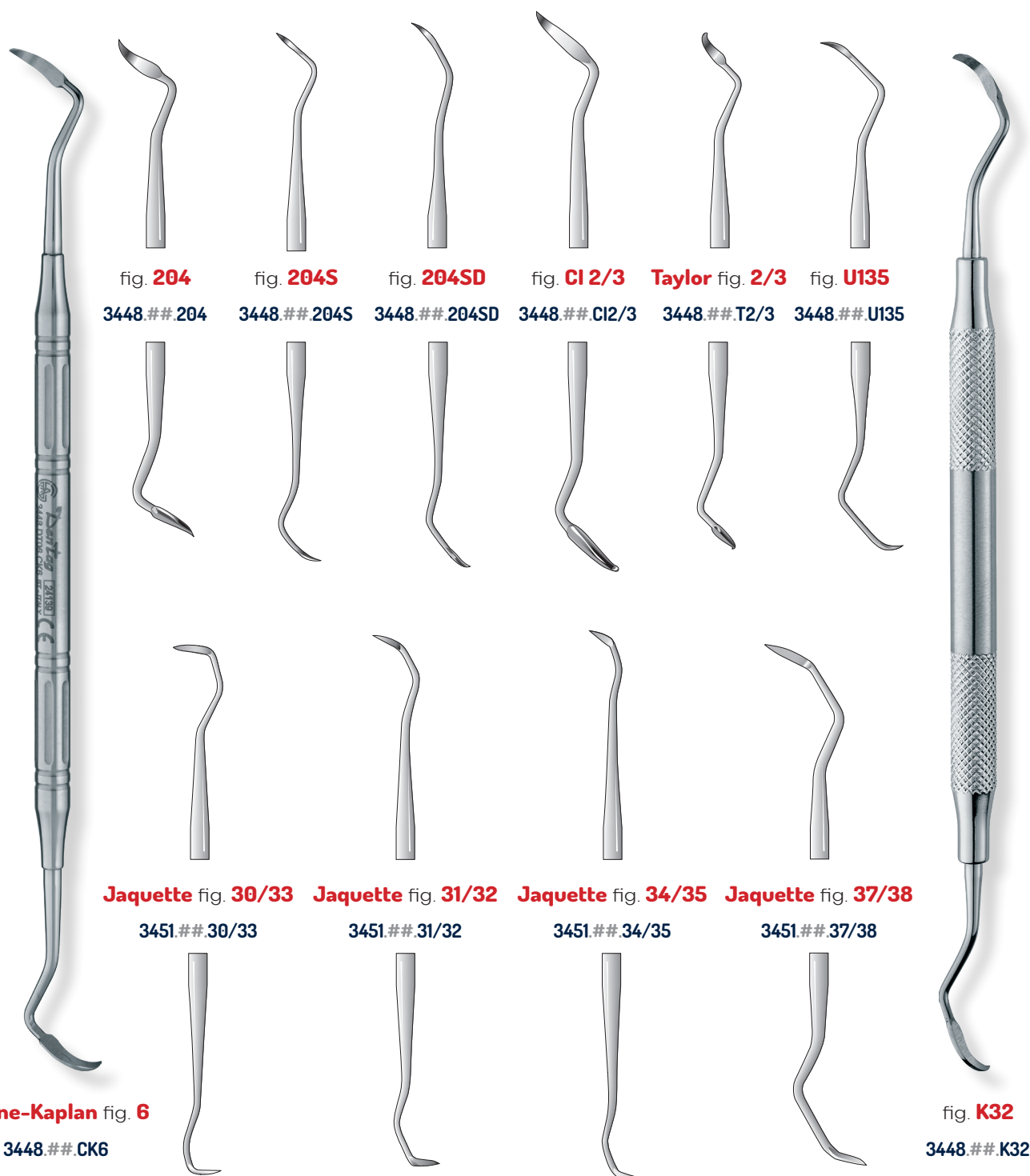
Available tips in  
Stainless Steel • Titanium • DLC

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